


9999 Businesspark
Ave, Suite B
San Diego, CA 92131
1-888-PSYONIC (779-
6642)
sales@psyonic.io


Socket Setup Guide (EQD)

User Manual: available on website for download at www.psyonic.io >
[Support](#) > [Clinicians](#)

Power Switch

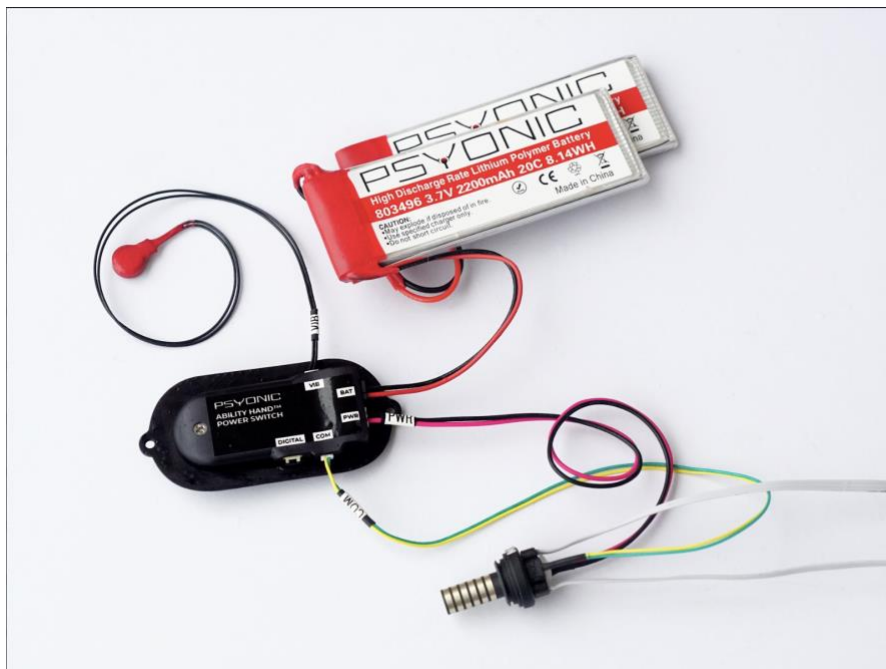

PSYONIC

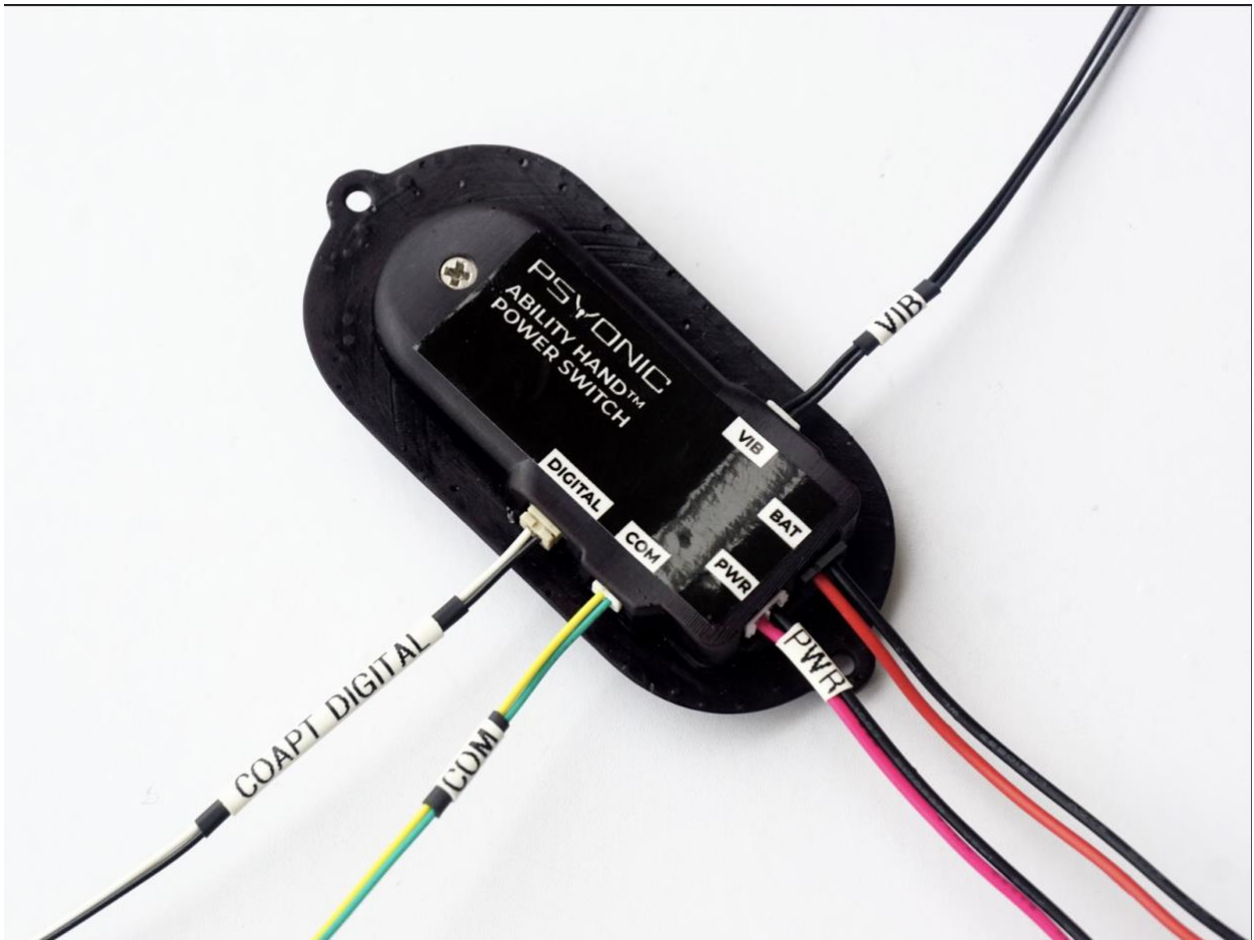


- The charging port should be situated proximally on the socket with the PSYONIC logo being readable from the user's viewpoint
- On the socket, the power switch should be positioned flat on surface (as much as possible) with the PSYONIC logo readable from the user's perspective
- Of the four (4) ports on the bottom of the power switch, the sticker provides a guide for cable connections. As a key:
 - VIB = PSYONIC Vibration Unit
 - BAT = PSYONIC Battery Pack
 - PWR = red/black cable; delivered already plugged in (POWER OUT)
 - COM = green/yellow cable; delivered already plugged in (Communication)
 - DIGITAL = black/white cable; delivered *only* for Coapt users
- The Power Switch can be positioned anywhere on the socket, per the prosthetist/technician design, though distal on a medial wall is recommended, so that the light on the switch is viewable while in use.

**Note: COAPT DIGITAL cable is only necessary for Coapt users.*

PSYONIC


Coaxial Plug

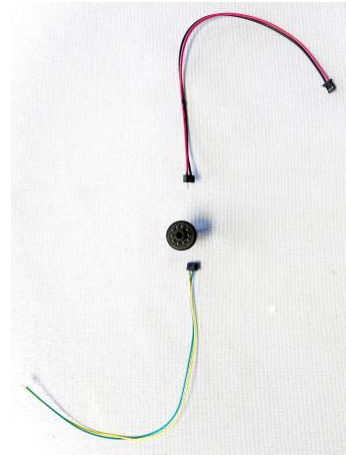

- PWR OUT (Power Out): red/black cable from Power Switch to be plugged into the 3-pin half moon adapter into the "1" port on coaxial plug



- - COM (Communication): green/yellow cable from Power Switch to be plugged into the 2-pin half moon end into the "A" port on coaxial plug

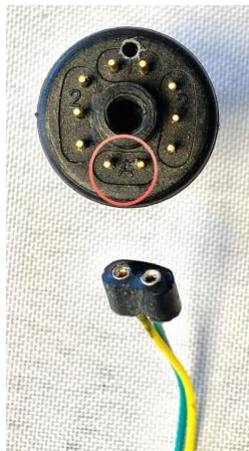

PSYONIC

Vibration Unit

- Optimal placement within the socket, between the inner socket and outer shell
- Not contacting user's skin directly

Battery Pack

- Packaged as a set of two 3.7V cells in series (7.4V, 2200mA total)
- Need to be plugged-in per the instructions under "Power Switch" section
- We recommend placement on a flat surface with the use of epoxy to adhere the battery pack within the inside of the socket

Charging

- Estimated time to charge: 1 hour to 100% with provided USB-C charger
- Cable can be inserted in any polarity direction
- Ports are USB-C and compatible with other USB-C cables; however, only the charging cable provided by PSYONIC has been tested to comply with the device for fast charge capability

Provided by PSYONIC:

- The Ability Hand
- 6-pin coaxial plug
- PSYONIC Battery Pack
- PSYONIC Power Switch
- USB-C 30W Wall Charger
- USB-C Charging Cable

What needs to be purchased separately:

- For the electronic quick disconnect (EQD) version of the Ability Hand
 - Coupling Piece
 - Ottobock 10S4 or Steeper QDACP
 - Lamination Ring
 - Ottobock 10S1=40/45/50/54 or Steeper QDALR 40/45/50/54
- For Myoelectric Direct Control
 - Electrodes
 - Ottobock, Steeper, or Motion Control
 - Electrode Cables
- For Pattern Recognition
 - Coapt Gen 2 System